www.photogensystem.com

www.photogensystem.com
www.mesoestetic.com

Grupo Body_esthetic Laboratories, S.L.
Av. Bertrán i Güell 72
08850 - Gavà (Barcelona)
Tel. + 34 902 26 20 31


made in Spain

mesoestetic® | medical devices

2953730903

Photogen System

Synergie Parfaite
Perfect Synergy

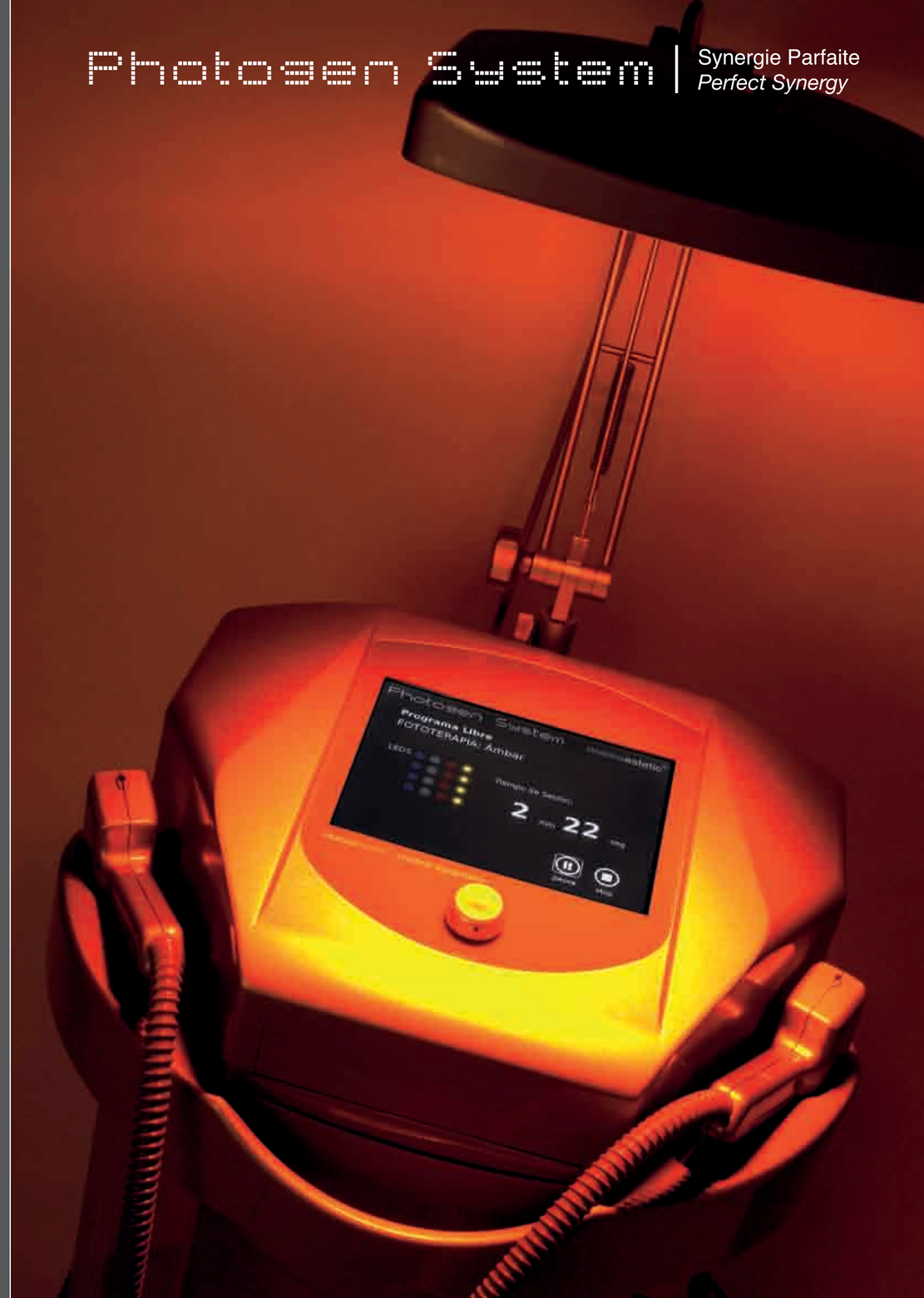Photosem
System

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Photogen System: is all benefits

- ***Smart combination of technologies:*** *Unipolar RF + Unipolar Transfer + Bipolar RF + Phototherapy*
- ***Technology Synergy + Medical Cosmetics***
- ***TriPhase** treatment **Protocols** to achieve complete treatments with spectacular results through user friendly programs*
- *Free and pre set **Programs***
- ***Multiple** Treatment Application*
- *Maximum **comfort, safety** and treatment **efficacy***
- *Long lasting and quick **results***
- ***Lifting effect without surgery***
- ***No consumables***

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Technical specifications

RF

Frequency	1 MHz
Energy density	150 J/cm ^{2}
Spot diameter	3 cm

PHOTOTHERAPY

Diodes number	140
LED Power	0,5 W/LED

Photogen System | RadioPhotonic Multitherapy by **meso**estetic®

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
Multiple Treatment Platform, Optimum results


***Photogen System (Radiophotonic Multitherapy)** synergistically combines the most effective Technologies **(Radiofrequency and Phototherapy)** with **Medical Cosmetics** formulated to increase the effects achieving spectacular results from the very first session.*

*Photogen System is a versatile medical-beauty platform with multiple body and facial treatments. It has **4 treatment programs** that are designed to meet the specific needs of each individual patient. It also includes a **free** mode program that enables the physician to customize treatments, which gives maximum **flexibility** for the patient.*

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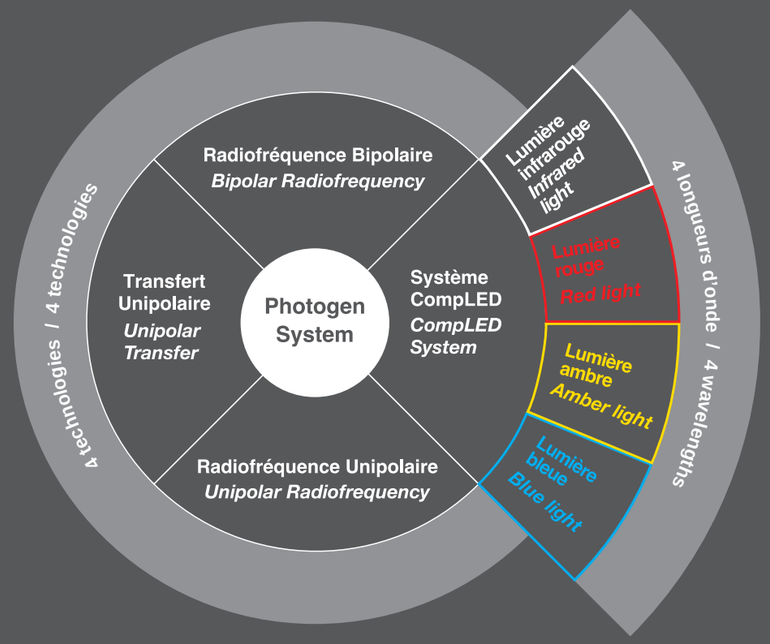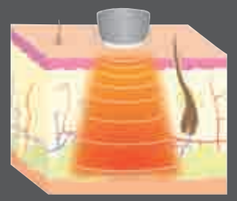
1. Unipolar Radiofrequency

When applied over the skin, RF energy produces a deep and selective heating into the dermis, which produces the contraction of collagen causing an immediate tightening effect.

*Over time, **dermal restructuring** is achieved. Substituting old collagen and elastin fibres with new ones, creating more elastic tissues, thus improving the skin's quality and extending the lifting effect.*

Unipolar radiofrequency acts locally in the deep dermal tissues, respecting the dermis and targeting skin flaccidity safely and effectively.

It is a very effective technique for the treatment of cellulite, localized fat, and intense body tightening.

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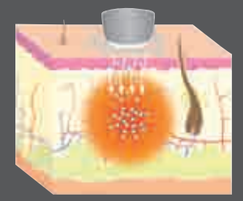2. Unipolar Transfer

*Photogen System introduces a revolutionary technology for the application of functional active principles. The Unipolar transfer applicator produces changes in the permeability of the skin due to its frequency settings, allowing **active principles to penetrate** into the dermis, where the repair cell mechanisms are found.*

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Immediate results.
Painless and non-ablative treatment.

Photogen System achieves spectacular results, painless, non-ablative and with no downtime. The patient can resume normal activity immediately post treatment.

The synergy achieved by the smart combination of the most innovative and efficient techniques, along with the use of revolutionary cosmetic formulations, brings out visible and long lasting results from the first session.

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
4 technologies

For the first time, comes a medical device combining 4 innovative and effective Technologies that work cohesively together, utilizing different levels to achieve optimum results:

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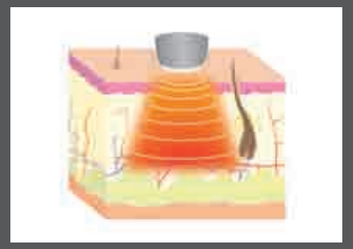
3. Bipolar Radiofrequency

RF is transmitted from an emitting pole to the receptor pole, which is found on the same handle. It acts more superficially than Unipolar Radiofrequency (0,5-1 cm), as its effects are limited to the area comprised within the two handheld poles.

RF energy produces a controlled increase of the dermal temperature, increasing the activity of the fibroblasts, which multiplies the collagen and elastin production, achieving an intense lifting effect without invasive surgery.

The treatment is completely safe, as the radiofrequency energy is applied selectively over the dermis avoiding thermal injury to the epidermis.

Optimum results on facial treatments and stretch marks treatments.

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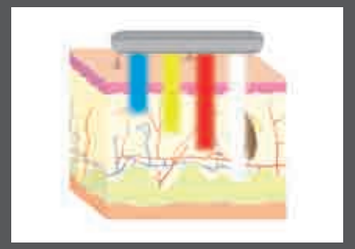4. High powered LED phototherapy: CompLED System

*Photogen System incorporates the most **complete** phototherapy treatment, allowing the penetration of electromagnetic energy through different wavelengths into the dermal layers.*

Cell structures absorb this energy, which have different benefits through the different deep levels.

*Its principles are based on **photodynamic cell stimulation** through visible Light produced by LEDs, which **inhibits or stimulates the different cell processes.***

*It is a **non-invasive** technique, which activates the cells through a process called "photo-modulation", with no contraindications.*

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HOME

For long lasting results of the cabin treatment, Photogen treatment packs include a maintenance cream for home use specially formulated with specific and high concentration of active ingredients.

Photogen System enables the therapist to work on free mode, adapting the working settings depending on the specific needs of each patient.

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10 sessions

In just 10 treatment sessions, Photogen System achieves spectacular results.

Treatment cycle
- *First 4 sessions (one session per week)*
- *Following 4 sessions*
 (one session every two weeks)
- *2 maintenance sessions (once per month)*

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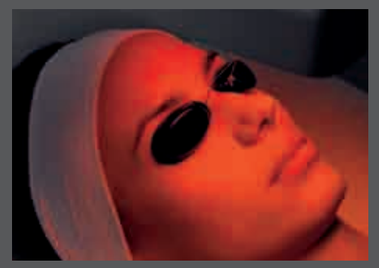
Application mode

Triphase preset protocols are easy to use and very effective from the first session, due to its diverse and unique multi-application treatment.

PHASE 1

*Active Principles Penetration on the dermis through the **Unipolar Transfer RF.***

PHASE 2

*Collagen and elastin production stimulation by using either **Unipolar or Bipolar RF,** depending on the depth of the area to be treated. To ensure optimal energy transmission, **meso**estetic® laboratory has developed the ADN Thermo Conducting gel extra moisturizing, which has hydrating and refreshing properties.*

PHASE 3

***High Powered Phototherapy.** Completes the treatment by means of specific wavelength LED light (red, blue, amber and infrared LEDs).*

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4 wavelengths

*The **CompLED system** allows the use of different wavelengths simultaneously which confers the unit to have maximum **versatility,** allowing different treatment combinations:*

Infrared light (880 nm):

*Intense regeneration, very effective for scar treatments. As it acts on deep layers, it is also very effective on cellulite and localized fat treatments. Acts on **deep dermal layers.***

Red light (660 nm):

*Produces a bio stimulation and deep repairing, stimulating the metabolism of the cells, multiplying fibroblasts activity during the production of elastin and collagen, and thus improving the quality of the tissues. Acts on deep **subcutaneous** tissues.*

Amber light (590 nm):

*Anti-aging and deep stretch marks treatments, increasing the luminosity of the skin and unifying skin tone. Also very effective on pigmented skin. Acts on **medium** tissues.*

Blue light (470 nm):

*Acne treatments. Acne bacteria are formed by a chemical substance called Porphyrine. The application of Blue Light over Porphyrine produces a chemical reaction, generating free radicals that destroy acne bacteria. Acts on **superficial** tissues.*

Photogen System includes different programs that allow the application of each wavelength individually or the different possible combination to customize treatments depending on the patient's individual needs.

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4 exclusive treatment programs

***meso**estetic® laboratories, leaders in the medical cosmetic field, have formulated four exclusive treatment programs to meet the demanding needs of each patient.*

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hydrating and anti wrinkle

Hydrophotonic solution program aims at hydrating the skin deeply, decreasing the depth of expression lines caused by hydration loss on the stratum corneum. It prevents the aging process by maintaining optimal skin hydration.

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facial tightening

Tightening Skin Solution program is specifically designed to tighten skin by restructuring the conjunctive tissue due to the remodelling of the collagen and elastin fibres, achieving a long lasting lifting effect without invasive surgery.

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cellulite and localized fat

Bi-Sculpting Solution program is specifically designed to remodel the shape of the body by helping reduce the appearance of localized fat, cellulite and skin flaccidity.

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body tightening

Bi-Tonic Solution program is specifically designed to reduce the appearance of body flaccidity through the restructuring of conjunctive tissue, due to the formation of new collagen and the contraction of collagen and elastin fibres.